



## TAURI7 Set Up Guideline



The Tauri Temperature-Check Tablets are designed to detect temperature anomalies. It's important to understand there are many factors, including environmental and physiological that can impact a person's surface temperature reading. Skin surface temperature vs actual core body temperature may differ either way. The Tauri Temperature-Check Tablets must be operated in accordance with the manufacturer's user guide. Tauri Temperature-Check Tablets are not intended nor designed to diagnose or detect medical conditions including, but not limited to, viruses or other illnesses. The Tauri temperature-check tablets should only be used to detect variations in surface temperature. In the event that an elevated skin temperature is detected, the subject should be advised to check temperature with an approved medical thermometer and the finding confirmed. Absence of an elevated skin temperature does not preclude a fever.

# Content

1. Initial Setup-----	2
2. Tools Included-----	3
3. Measuring Distance -----	4
4. Mounting Instructions-----	5-13
5. Connecting to Wi-Fi -----	14-15
6. Setting the Time Zone-----	16 -17
7. TAURI Support-----	18

## 1. Initial Setup

Connect the cable from power adaptor to the DC in of TAURI 7

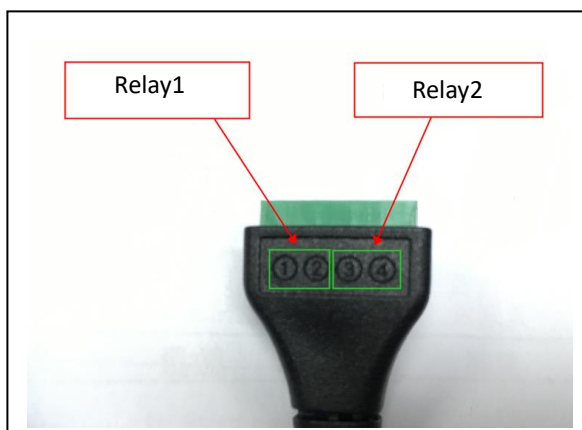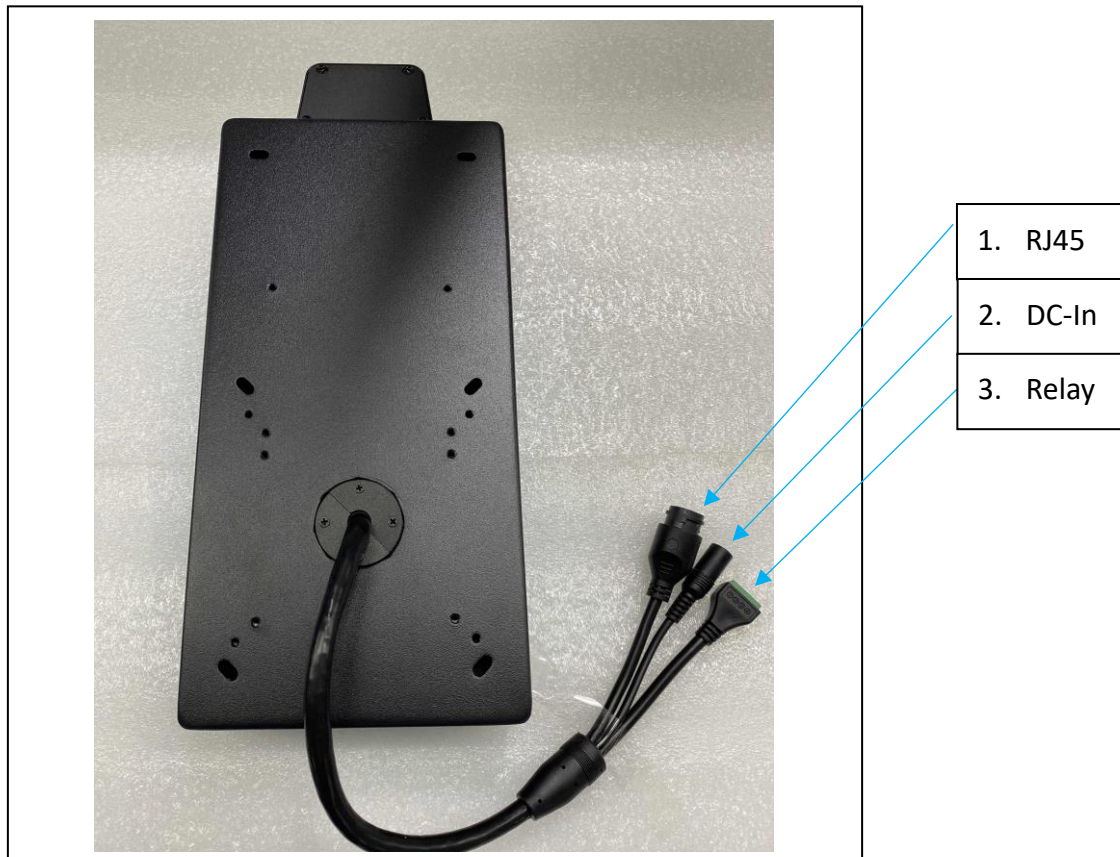


Relays




Power adaptor

## 2. Tools Included

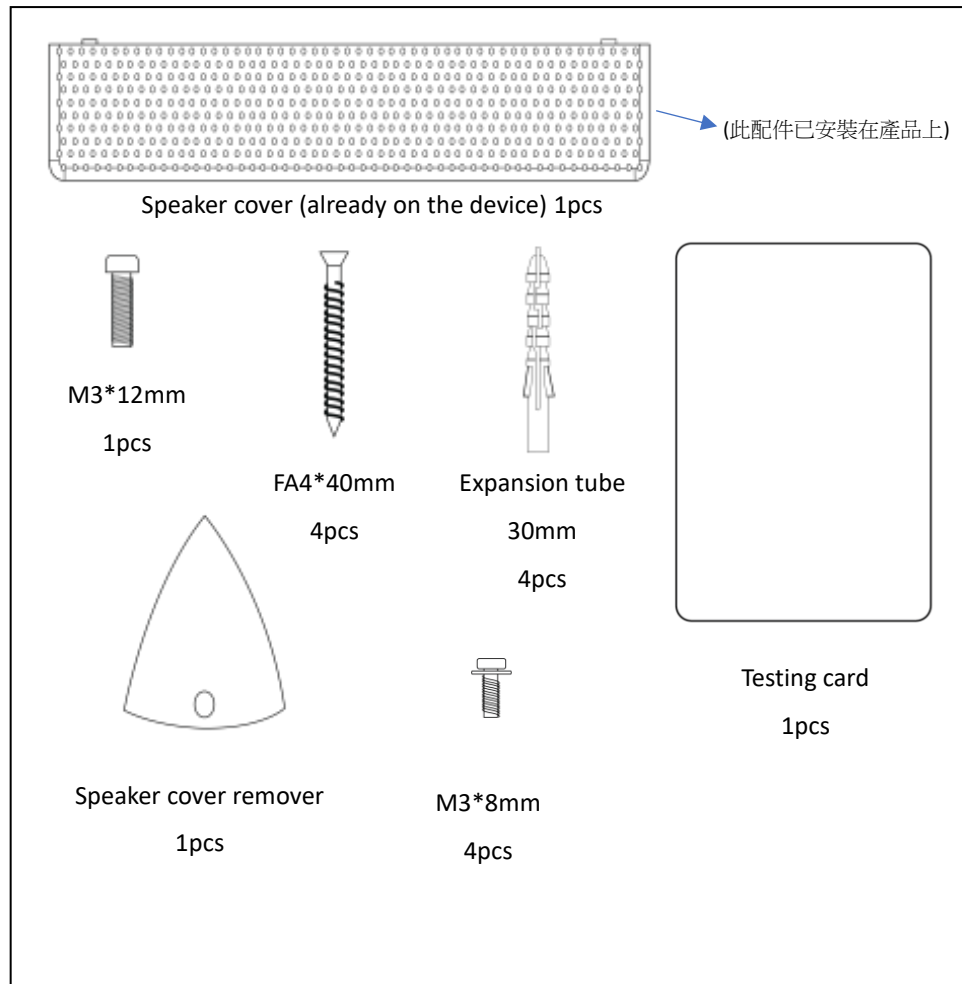

### 3. Measuring distance

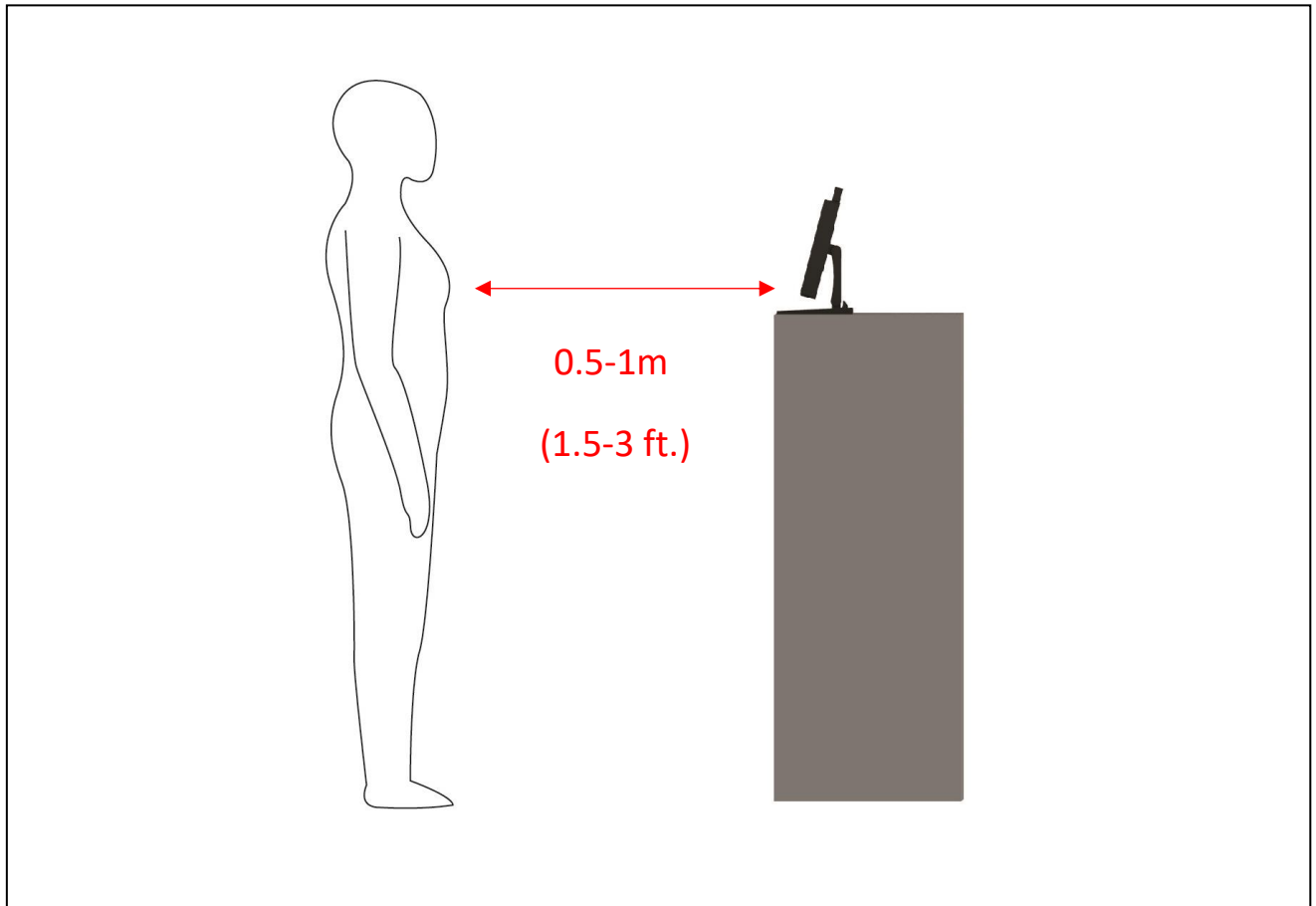

1. Please stand within 0.5 to 1m (1.5-3 ft.) from TAURI7 for best result.

## 4. Mounting Instructions


### Installation Example – TAURI7 with FIX028 table stand




1. Small size mount.




2. Move the arm in vertical direction.

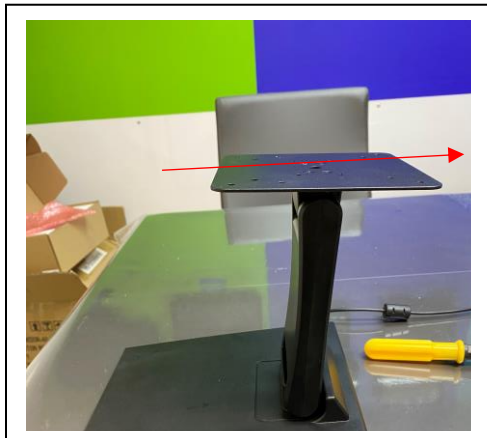


3. Adjust the plate in horizontal position.

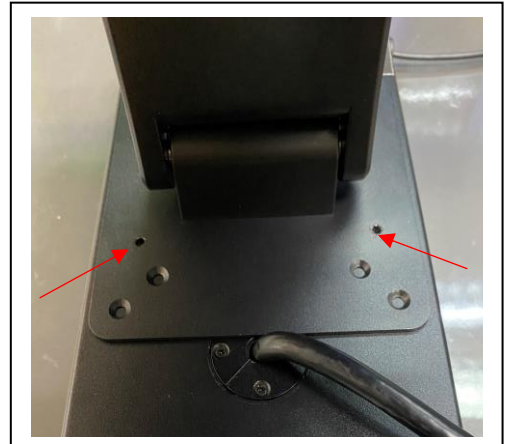

4. Place TAURI7 on table and screw two 3M x 8mm screws here.

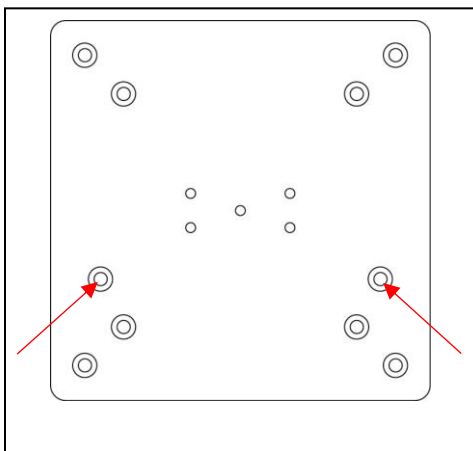



5. The position of the holes in step 4 on metal plate.

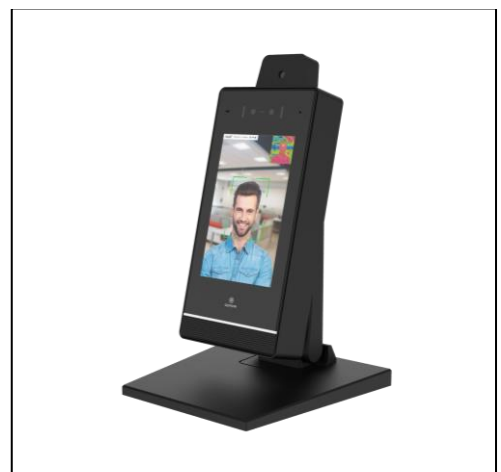

6. Finish.

## Installation Example – Installing TAURI7 using wall mount


**\*Please note: When mounting TAURI7 to wall, PoE is the only recommended method for power.**

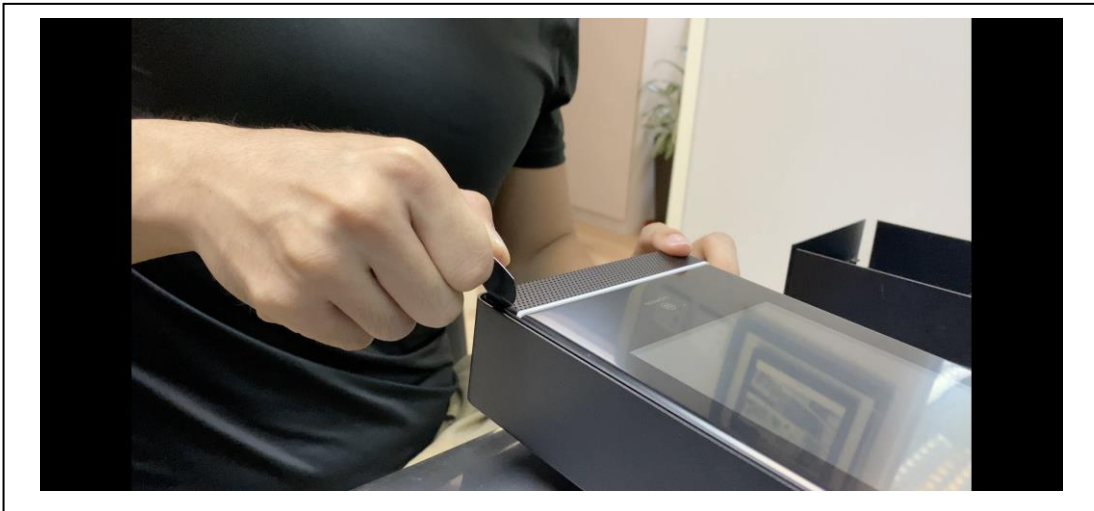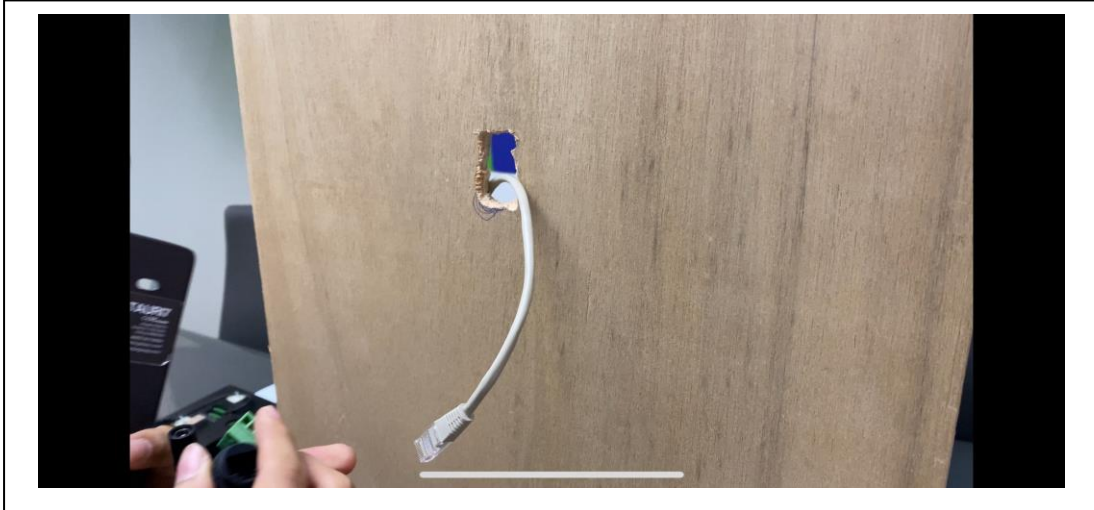


1. Remove speaker cover by using speaker remover.

Place the speaker cover remover perpendicular to the edge of TAURI7.  
Gently lift the speaker cover up to remove it from TAURI7.



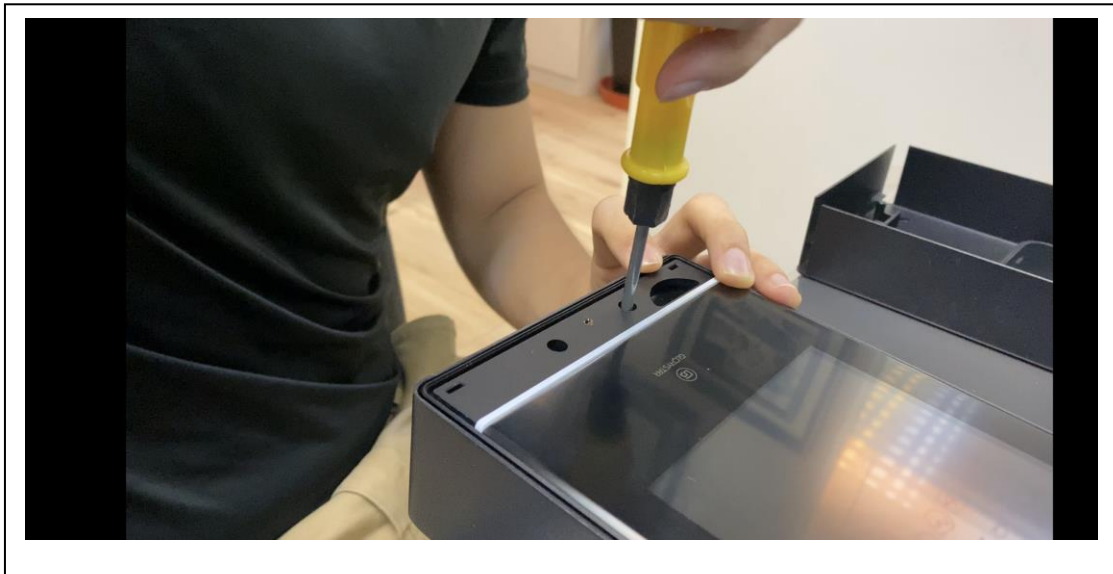


2. Release the screws under the speaker cover and take out device from the metal case.

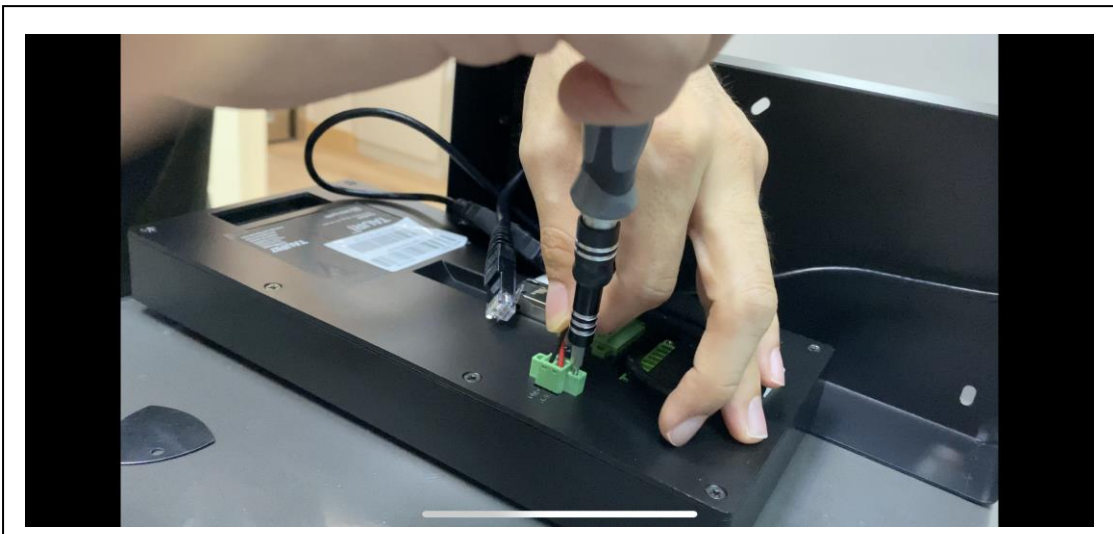


3. Release all cables and USB on the back.




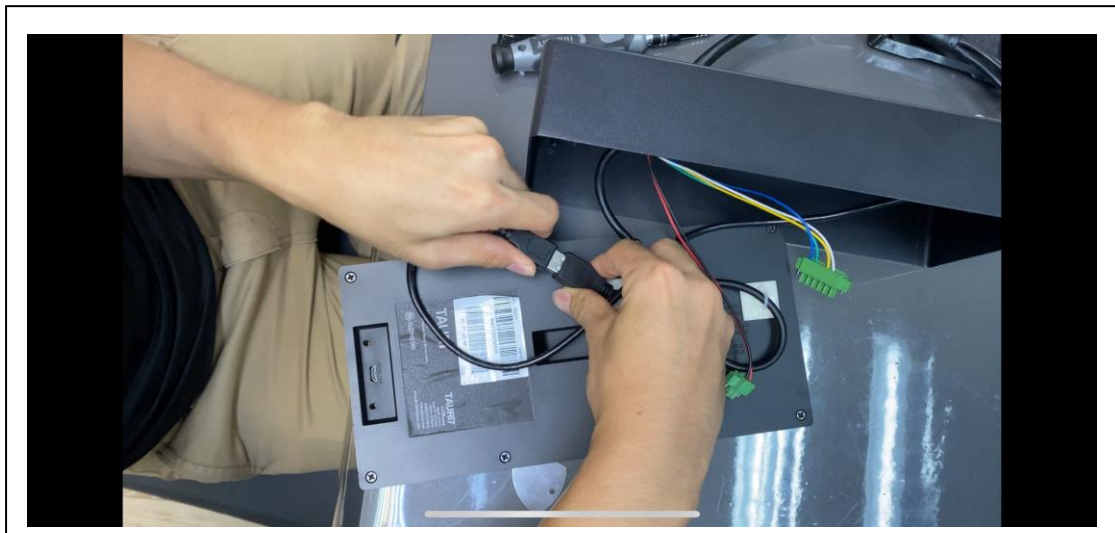


4. Separate device and metal case.

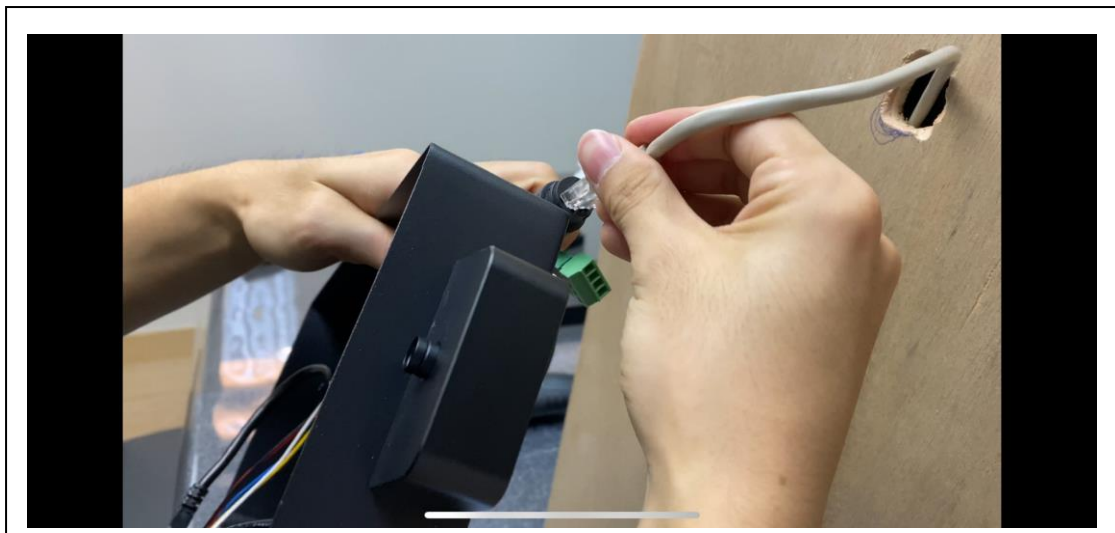

5. Connect PoE from wall to the PoE on metal case.

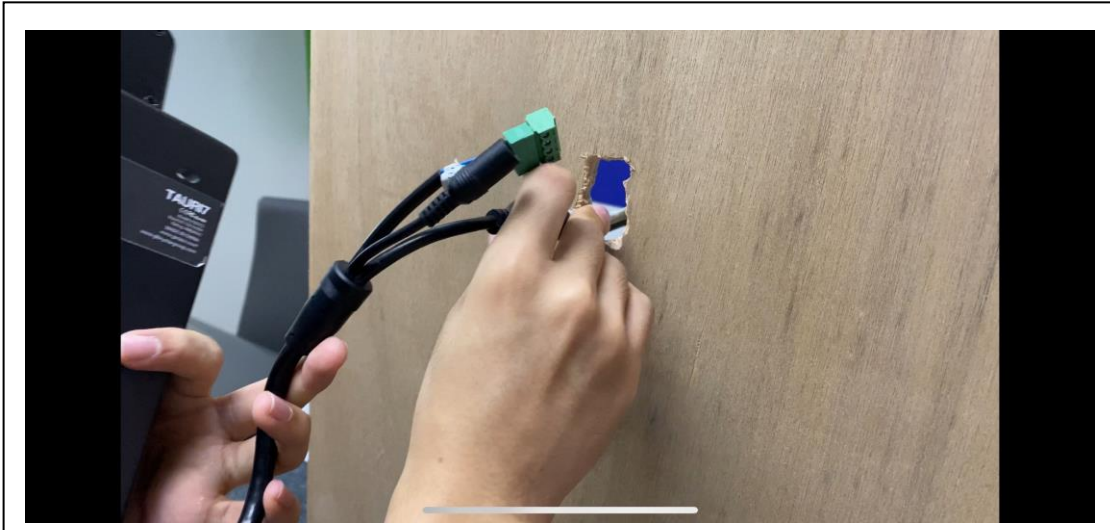


6. Put all cables into the wall.

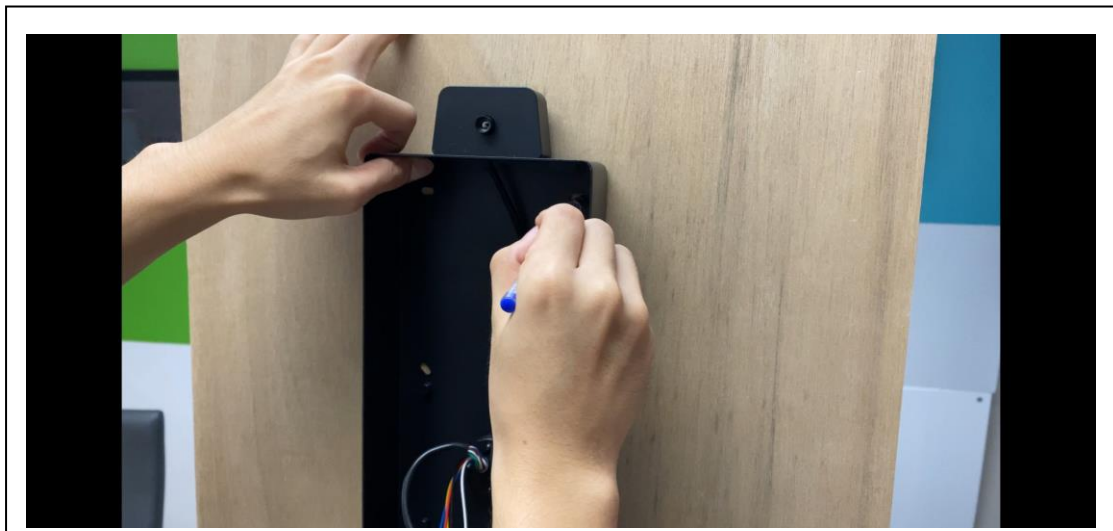


7. Place metal case on wall and mark the four corner hole positions.



6. Drill four holes on the wall.

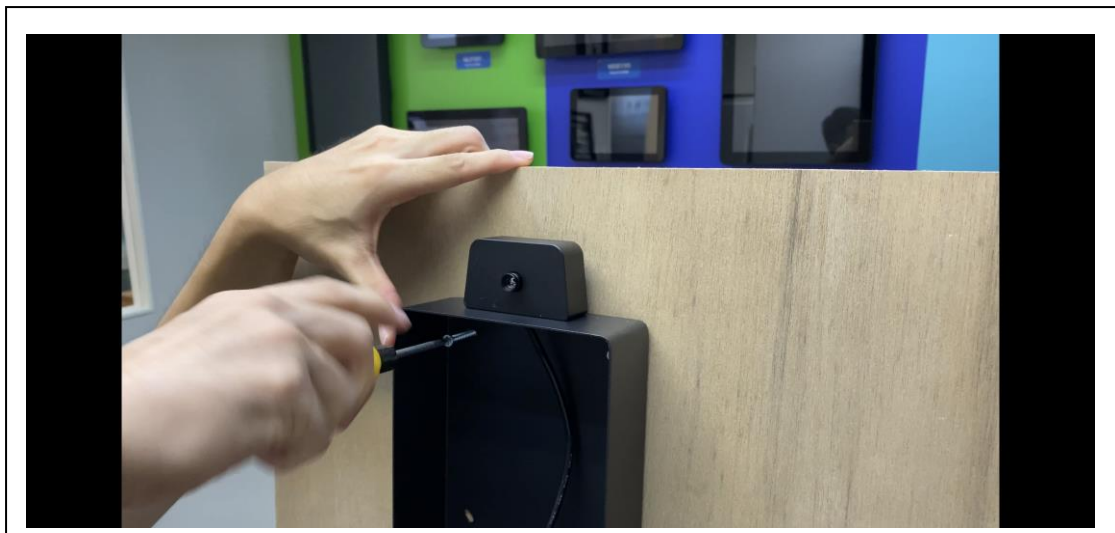


7. Screw 40mm screws on four holes on the metal case.  
(You may need to insert expansion tube for concrete wall.)

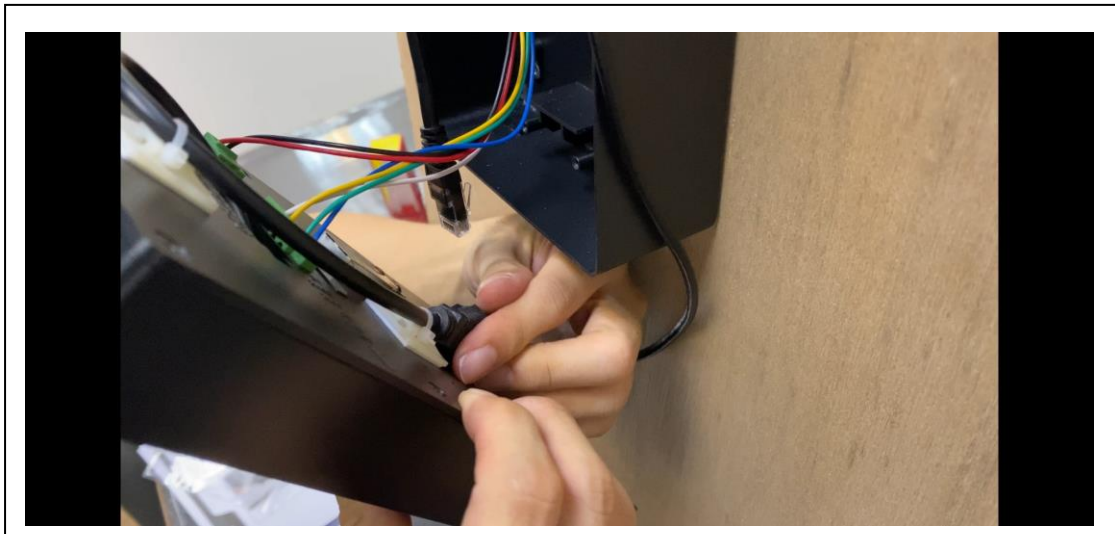

8. Connect cables and USB then insert device into the metal case.



9. Screw two front screws.

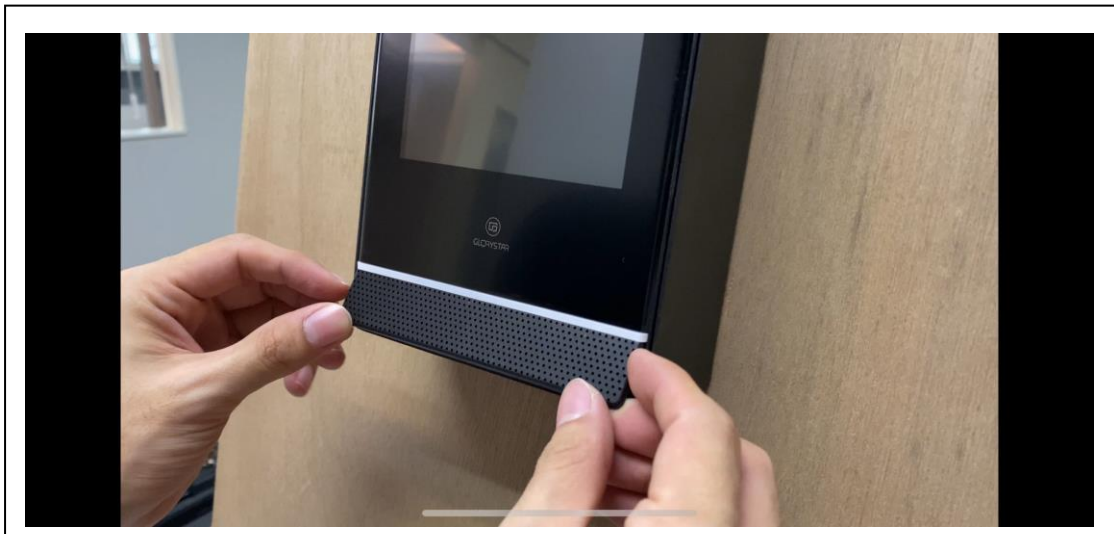

11. Put the speaker cover on.



12. Finish.




**Please scan the QR code for demo video**




## 5. Connecting to Wi-Fi

**Please connect Wi-fi before starting set-up**


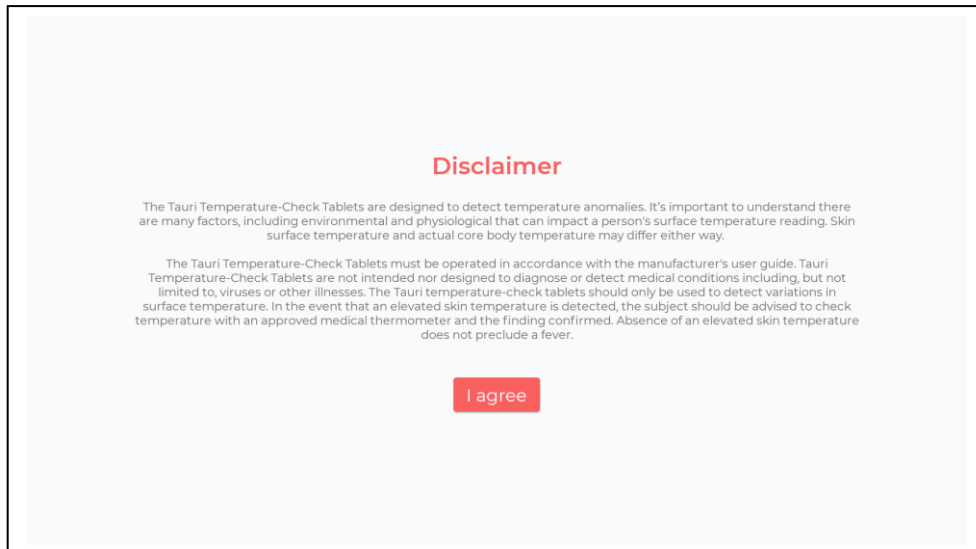


5.1 Press “I agree” and start software setup.

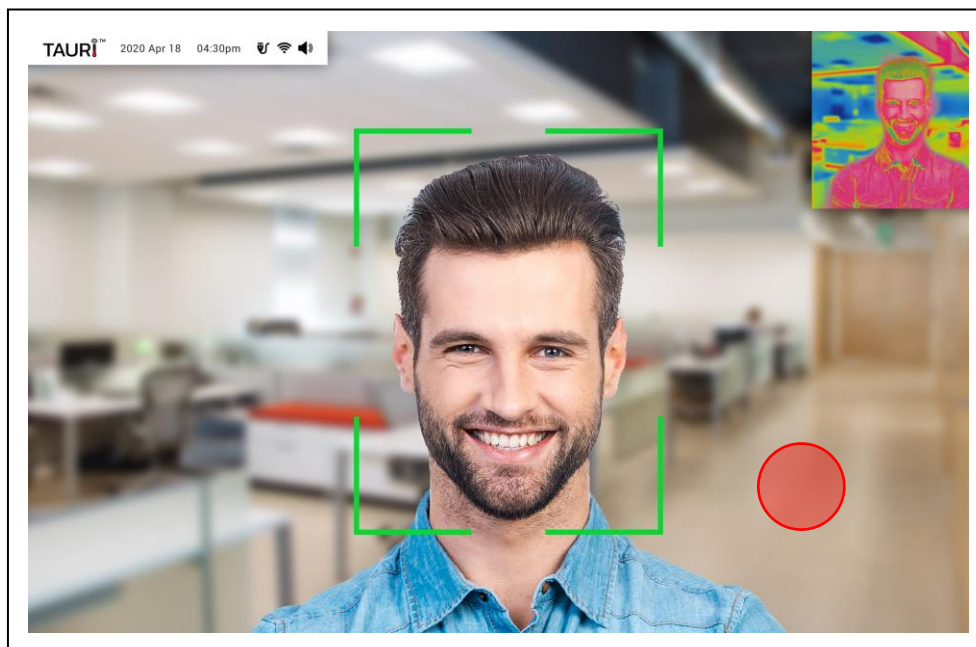


5.2 Press and hold on screen for 5 seconds and enter “0000” to access into settings. Then select “system settings” and connect Wi-fi.



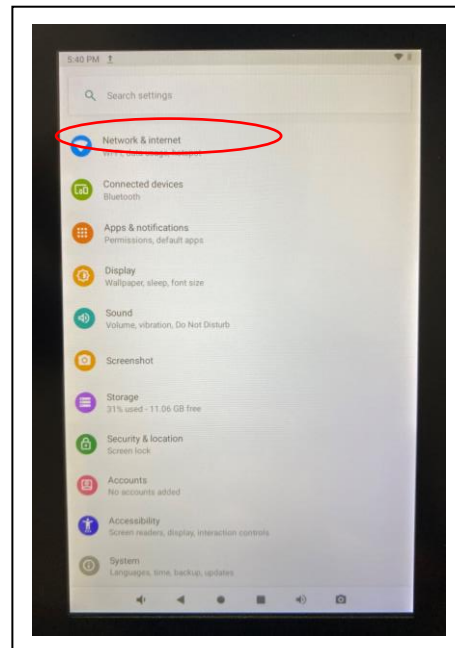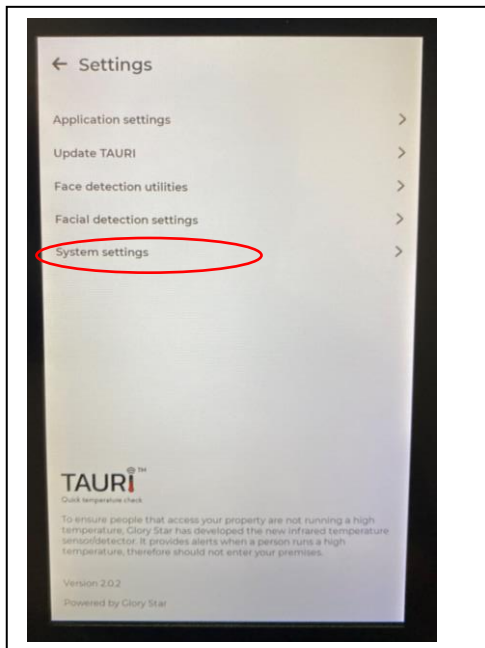


### 3.3 Select “system settings” and connect Wi-if

## 6. Setting the Time Zone

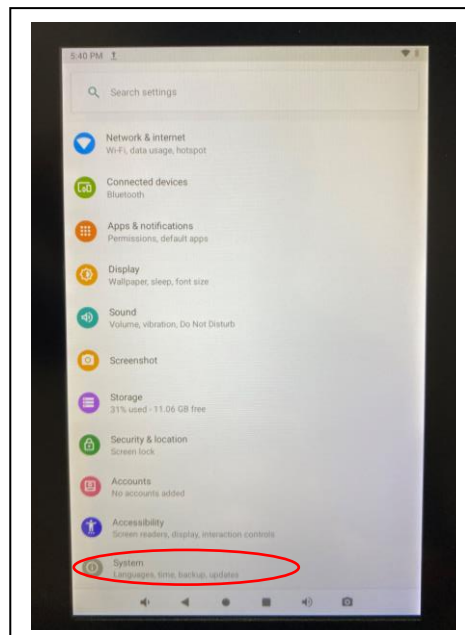


6.1 Select "System settings" in TAURI main page.

6.2 Select "System".

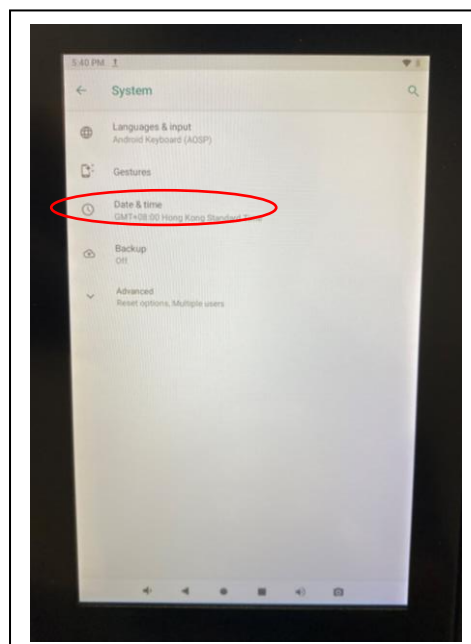

4.2. Select "Date & time"

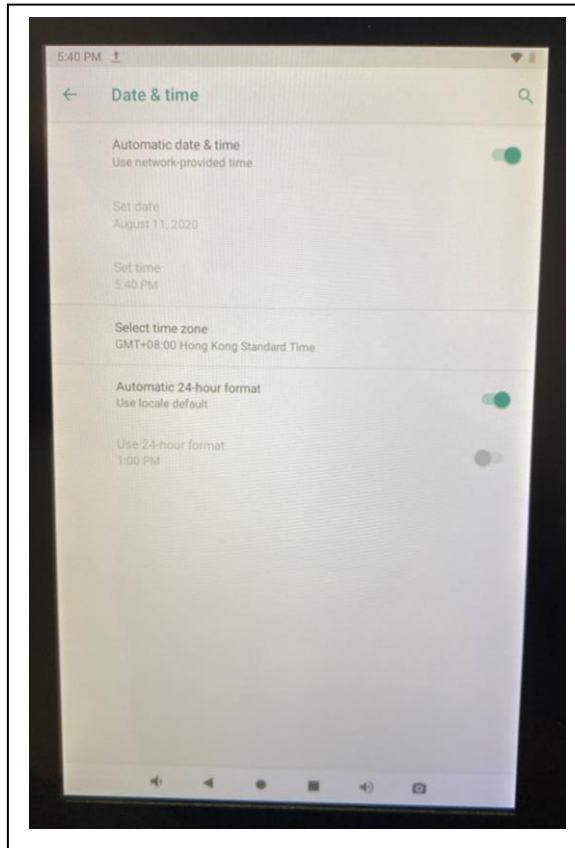

4.4. Select your time zone.

## **TAURI Support**

Flat D, 5/F, Excelsior Building, 364 Nathan Rd., Kowloon, H.K

Hong Kong: (852) 2376 2803

USA: (626) 228 5927

Website: [www.gotauri.com](http://www.gotauri.com)